




Incentive spirometry

SPIRO-BALL® incentive exerciser and TRI-BALL® respiratory exerciser

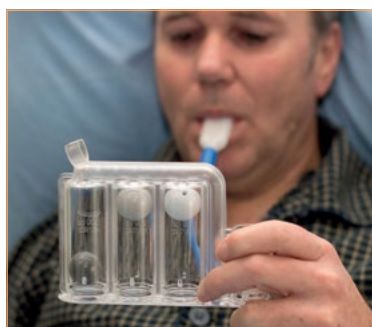

Oxygen and Aerosol Therapy • Respiratory Physiotherapy

SPIRO-BALL®

SPIRO-BALL® Incentive Exerciser is a single patient use device designed to help the patient to recover normal breathing with an indicator that allows monitoring progress.

Inspiration tube

Flexible tube with mouthpiece

Two particle filters in the inspiration system

Minimise the risk of inhaling foreign bodies

Flow indicator

Intuitive flow control with a smiling face

Target pointer

Promote the exercise with the target pointer



| Code | Description | Box Qty. |
|---------|--|----------|
| 1537000 | SPIRO-BALL® Incentive Exerciser 4000ml (adult) | 12 |

TRI-BALL®

TRI-BALL® encourages the patient to achieve a maximum prolonged inhalation, named SMI technique (Sustained Maximal Inspiration). The SMI is an exercise of deep, slow breathing that when repeated several times per day, can prevent respiratory diseases.

Flexible tube with mouthpiece

Includes a particle filter

Three coloured balls in independent chambers

Allow exercise of three different levels of inspiration

Single patient use

Reduces the risk of cross-contamination



| Code | Description | Box Qty. |
|---------|---------------------------------|----------|
| 1538000 | TRI-BALL® Respiratory Exerciser | 12 |

Please note these products are only available in select markets.

For further information, contact your local office.

SPIRO-BALL® and TRI-BALL® are registered trademarks of Leventon S.A.U.

Please think before you print Save energy and paper. If you must print this information sheet please print it double sided.

IS10.14 • Issue 4 01.23



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 9656 300 F: +44 (0)118 9656 356 info@intersurgical.com www.intersurgical.com

UK • Ireland • France • Germany • Spain • Portugal • Italy • Benelux • Sweden • Denmark • Lithuania • Russia • Czech Republic
 Turkey • South Africa • China • Japan • Taiwan • Philippines • USA • Canada • Colombia • Australia