




STARMed CaStar R Up

Openable hood for non-invasive mechanical ventilation (NIV)

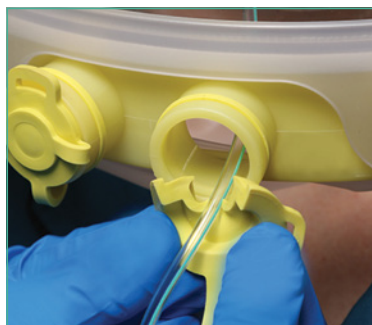

Critical Care ▪ Patient Interfaces ▪ **STARMed**



Quality, innovation and choice

Interact with us



www.intersurgical.co.uk



STARMed CaStar R Up openable hood for NIV

StarMed CaStar R Up

Openable CaStar R Up is a comfortable, versatile and lightweight patient interface for non-invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R Up is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in semi-recumbent, supine and prone positions providing an ideal option for NIV therapy.

Advantages

- Allows rapid access to the patient during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood

Technical Features

- Manufactured from phthalate-free materials
- 22M connectors compliant to ISO standard
- Built in inflatable neck cushion for added patient comfort
- Comfortable and atraumatic underarm fastening straps
- Other fastening accessories are available if required
- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Read more evidence on www.intersurgical.co.uk



CaStar R Up openable hood for NIV with bi-directional anti-asphyxiation valve and patient access port, inflatable neck cushion and inflation bulb, ear plugs and measuring tape*

Code	Size	Neck circumference	Box Qty.
CA022XXL/2RW	Extra Extra Large	≥50cm	5
CA022XL/2RW	Extra Large	45-52cm	5
CA022L/2RW	Large	40-47cm	5
CA022M/2RW	Medium	34-41cm	5
CA022S/2RW	Small	27-34cm	5
CA022XS/2RW	Extra Small	17-27cm	5


www.intersurgical.co.uk/info/starmed
Make an enquiry

*Tubing not included and has to be ordered separately. Please refer to the website for the range of accessories.


 Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: 0118 965 6300 F: 0118 965 6356 info@intersurgical.co.uk www.intersurgical.co.uk


 Management System
 ISO 13485:2016
 ISO 9001:2015

www.tuv.com
 ID 900005840

Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print
Save energy and paper. If you must print this information sheet please print it double sided.

IS1.39_UK • Issue 5 06.20