




STARMed CaStar Up

Openable hood for CPAP therapy

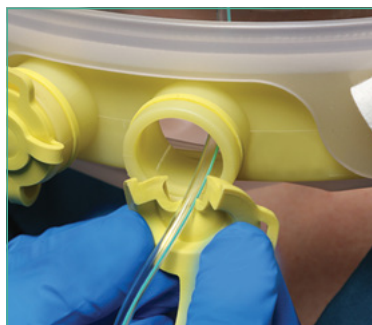

Critical Care ▪ Patient Interfaces ▪ **STARMed**



STARMed CaStar Up openable CPAP hood

StarMed CaStar Up

Openable CaStar Up is a comfortable, versatile and lightweight CPAP patient interface designed to offer an alternative to traditional methods of CPAP delivery. The CaStar Up is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for CPAP therapy.

Advantages

- Allows rapid access to the patient during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood

Technical Features

- Manufactured from phthalate-free materials
- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Built in inflatable neck cushion for added patient comfort
- Comfortable and atraumatic underarm fastening straps
- Other fastening accessories are available if required
- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Read more evidence on www.intersurgical.co.uk

Redesigned patient access port
With bi-directional anti-asphyxiation valve

Transparent
For ease of patient monitoring

Integrated pressure manometer

Easy removal
For rapid patient access

New design sealed access ports
Closer to the patient's airway

PHF
PHTHALATE-FREE

CaStar Up openable CPAP hood with bi-directional anti-asphyxiation valve and patient access port, inflatable neck cushion and inflation bulb, ear plugs and measuring tape*

Code	Size	Neck circumference	Box Qty.
CA012XXL/2W	Extra Extra Large	≥50cm	5
CA012XL/2W	Extra Large	45-52cm	5
CA012L/2W	Large	40-47cm	5
CA012M/2W	Medium	34-41cm	5
CA012S/2W	Small	27-34cm	5
CA012XS/2W	Extra Small	17-27cm	5



www.intersurgical.co.uk/info/starmed



Video available at www.intersurgical.co.uk

[Make an enquiry](#)

*PEEP valve and tubing not included and have to be ordered separately. Please refer to the website for the range of accessories.



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



Manufacturer: Intersurgical S.p.A. - Mirandola (MO) - Italy
The Intersurgical S.p.A. Organisation has a Quality Management System Certificate by KIWA CERMET ITALY S.p.A. according to ISO 9001:2015 and ISO 13485:2016

Please think before you print
Save energy and paper. If you must print this information sheet please print it double sided.

IS1.31_UK • Issue 6 02.19