




STARMed CaStar Up

Openable hood for CPAP therapy

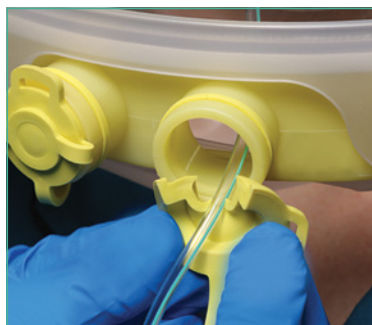

Critical Care ▪ Patient Interfaces ▪ STARMed



Quality, innovation and choice

Interact with us



www.intersurgical.co.uk



STARMed CaStar Up openable CPAP hood

StarMed CaStar Up

Openable CaStar Up is a comfortable, versatile and lightweight CPAP patient interface designed to offer an alternative to traditional methods of CPAP delivery. The CaStar Up is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in semi-recumbent, supine and prone positions providing an ideal option for CPAP therapy.

Advantages

- Allows rapid access to the patient during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood

Technical Features

- Manufactured from phthalate-free materials
- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Built in inflatable neck cushion for added patient comfort
- Comfortable and atraumatic underarm fastening straps
- Other fastening accessories are available if required
- Two versatile sealed access ports suitable for single or multiple lumen probes or catheters
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Read more evidence on www.intersurgical.co.uk

Redesigned patient access port
With bi-directional anti-asphyxiation valve

Transparent
For ease of patient monitoring

Integrated pressure manometer

Easy removal
For rapid patient access

New design sealed access ports
Closer to the patient's airway

CaStar Up openable CPAP hood with bi-directional anti-asphyxiation valve and patient access port, inflatable neck cushion and inflation bulb, ear plugs and measuring tape*

Code	Size	Neck circumference	Box Qty.
CA012XXL/2W	Extra Extra Large	≥50cm	5
CA012XL/2W	Extra Large	45-52cm	5
CA012L/2W	Large	40-47cm	5
CA012M/2W	Medium	34-41cm	5
CA012S/2W	Small	27-34cm	5
CA012XS/2W	Extra Small	17-27cm	5

www.intersurgical.co.uk/info/starmed
 Video available at www.intersurgical.co.uk

*PEEP valve and tubing not included and have to be ordered separately. Please refer to the website for the range of accessories.

INTERSURGICAL
 COMPLETE RESPIRATORY SYSTEMS

Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: 0118 965 6300 F: 0118 965 6356 info@intersurgical.co.uk www.intersurgical.co.uk

Management System
 ISO 13485:2016
 ISO 9001:2015

www.tuv.com
 ID 900005840

Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print
Save energy and paper. If you must print this information sheet please print it double sided.

IS1.31_UK • Issue 7 06.20