



# STARMed CaStar Ped

Paediatric hood for CPAP therapy and non-invasive mechanical ventilation (NIV)



Critical Care ▪ Patient Interfaces ▪ **STARMed**



Quality, innovation and choice

Interact with us



[www.intersurgical.co.uk](http://www.intersurgical.co.uk)



# STARMed CaStar Ped paediatric hood for CPAP therapy and NIV

## StarMed CaStar Ped

The CaStar Ped is a comfortable, versatile and lightweight CPAP and NIV patient interface for paediatrics. It is designed to offer an alternative to traditional methods of CPAP and NIV delivery. The CaStar Ped is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for CPAP therapy and NIV.


### Technical Features

- Manufactured from phthalate-free materials
- 22M/22F connectors compliant to ISO standard (CPAP version)
- 22M connectors compliant to ISO standard (NIV version)

- Fitted with integrated pressure gauge to check the presence of pressure within the hood (CPAP version only)
- Comfortable and atraumatic underarm fastening straps
- Two versatile sealed access ports suitable for probes or catheters of various gauges
- New design patient access port. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

#### Bibliography:

- Piastra M, et al. Noninvasive pressure-support ventilation in immunocompromised children with ARDS: a feasibility study, Intensive Care Med. 2009 Aug;35(8):1420-7
- Piastra M, et al. Treatment of acute respiratory failure by helmet-delivered non-invasive pressure support ventilation in children with acute leukemia: a pilot study, Intensive Care Med. 2004 Mar;30(3):472-6

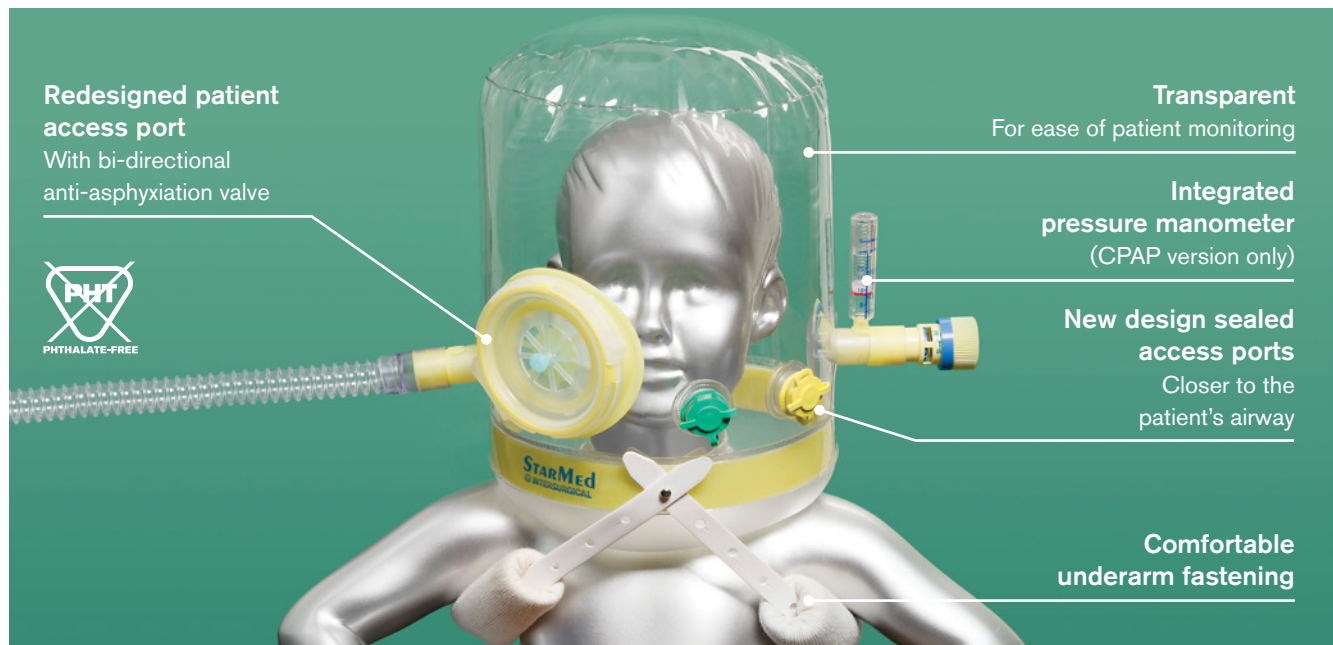

## CaStar Ped paediatric hood for CPAP therapy and NIV with bi-directional anti-asphyxiation valve and patient access port\*

Code	Description	Patient weight	Box Qty.
CP201PED/2R	CaStar R Ped for NIV with 22M/22M connections	>15Kgs	5
CP200PED/2	CaStar Ped for CPAP with 22M/22F connections	>12Kgs	5



[www.intersurgical.co.uk/info/starmed](http://www.intersurgical.co.uk/info/starmed)

Make an enquiry

\*PEEP valve and tubing not included and have to be ordered separately. Please refer to the website for the range of accessories.



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



Manufacturer: Intersurgical S.p.A. - Mirandola (MO) - Italy  
The Intersurgical S.p.A. Organisation has a Quality Management System Certificate by KIWA CERMET ITALY S.p.A. according to ISO 9001:2015 and ISO 13485:2016

**Please think before you print**  
Save energy and paper. If you must print this information sheet please print it double sided.

IS1.19\_UK • Issue 7 02.19